













## 50/15 PRIME E 10 BAR

**Cod. PR-0F200-040-AN**


Dry compressor, 2 cylinders, flow rate 110 l/min @ 8 bar, 40 lt tank equipped with an absorption dryer capable of delivering compressed air with filtration down to 0.01 micron and dew point below -40 ° C at atmospheric pressure.

 Maximum simultaneous use: 2

 Dental Chairs: 3

									
Air Dryer	FAD @ 8bar	Motor	Voltage Hz	Receiver	Noise Level	Duty Cycle	max P.	Size	net Weight
Yes	m <sup>3</sup> /min 0.11    l/min 110	kW 2.2	230V 50Hz	l 40	db(A) 66	100 %	bar 10	mm 740x410x770	Kg 47

## Air delivery performance

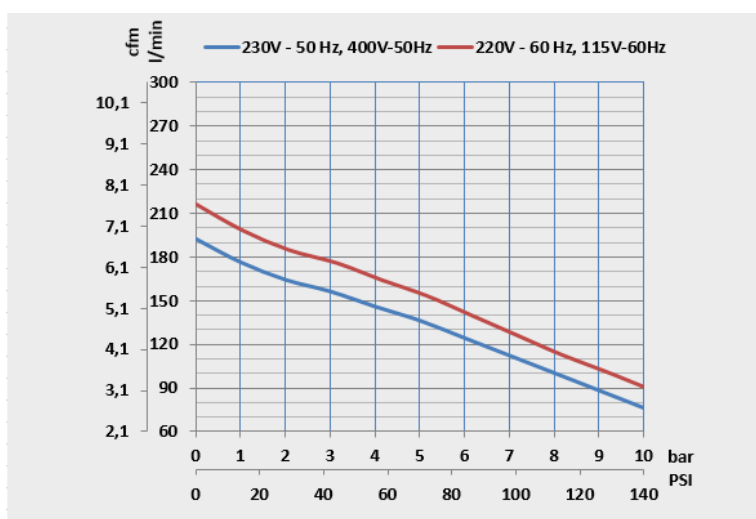

## Maintenance kits

Kit.	Description	Type	Hours
TA-FL-026-K-SCP	<p><b>PRIME 2-3 M filter cartridges kit</b></p> <p>Cartridge replacement kit for ordinary maintenance of the PURE AIR PRIME 2.3 M range.</p> <p>Comprehends:</p> <ul style="list-style-type: none"> <li>• 3 suction filters, code: TA-FL-058-C</li> <li>• 1 5 micron filter cartridge, code: TA-FL-001-C-SM</li> <li>• 1 0.01 micron filter cartridge, code: TA-FL-024-C2-SM</li> <li>• 2 o-rings for filters, cod. TA-FL-001-OR</li> <li>• 1 non-return valve pads, code: CO-PV-001</li> <li>• 1 5 micron filter cartridge (AF30), code: TA-FL-075-C</li> <li>• 1 0.01 micron filter cartridge (AFD30), code: TA-FL-076-C</li> <li>• 2 o-rings for filters, cod. TA-FL-075-OR</li> <li>• 1 AMG150 water separator filter cartridge, code: TA-FL-074-C</li> </ul> <p>Replace the filtration cartridges every 2 years.</p>	Ordinary	1.500h/2y
PR-GC-031-K-SCP	<p><b>PRIME 2</b></p> <p>Extraordinary maintenance spare parts kit PRIME 1.5 kW (PRIME 2) versions.</p> <p>comprehends:</p> <ul style="list-style-type: none"> <li>• 2 suction filters, code: TA-FL-058-C</li> <li>• 4 O-rings, code: MA-PC018</li> <li>• 2 piston rings , code: MA-MG00001</li> <li>• 2 piston guide rings, code: MA-MG00003</li> <li>• 2 expansion rings for the segments, code: MA-MG00024</li> <li>• 2 complete valve plates, code: MA-MG-02002P-C-P</li> <li>• 1 non-return valve pads, code: CO-PV-001</li> </ul> <p>The extraordinary maintenance intervals recommended for this compressor are 4000 hours/4 years.</p>	Extraordinary	4.000h/4y