



## CS 30/15 PRIME M

**Cod. PR-CS200-040-M**

Oil free reciprocating air compressors with soundproof cabinet, 2 cylinders, with a compressed air FAD of 170 l/min @ 5 bar pressure. Receiver capacity 40 L, with internal surface treatment to protect against corrosion. Complete with membrane dryer to ensure maximum air purity. This compressor is available also in **10 BAR VERSION**, visit the appropriate section!


 Maximum simultaneous use: 2

 Dental Chairs: 3

Medical Device CE 0425

Cod.	Air Dryer	Air capacity	FAD FAD @ 5bar	Motor	Voltage Hz	Receiver	Noise Level	RPM	Duty Cycle	max P.	Size	net Weig
CS 30/15 PRIME M PR-CS200-040-M	Yes	m³/min/l/min 0.23 230	m³/min l/min 0.17 170	kW 1.5	230V 50Hz	l 40	db(A) 53	1360	70 %	bar 7	mm 490x720x890	Kg 94
CS 30/15 PRIME M PR-CS200-040-M-400	Yes	m³/min/l/min 0.23 230	m³/min l/min 0.17 170	kW 1.5	400V 50Hz	l 40	db(A) 53	1470	70 %	bar 7	mm 490x720x890	Kg 94
CS 30/15 PRIME M 220V - 60Hz PR-CS200-040-M-6-UL	Yes	m³/min/l/min 0.28 280	m³/min l/min 0.185 185	kW 1.5	220V 60Hz	l 40	db(A) 53	1760	70 %	bar 7	mm 490x720x890	Kg 94
CS 30/15 PRIME M 110V - 60Hz PR-CS200-040-M-110	Yes	m³/min/l/min 0.28 280	m³/min l/min 0.185 185	kW 1.5	110V 60Hz	l 40	db(A) 53	1760	70 %	bar 7	mm 490x720x890	Kg 94

## Air delivery performance

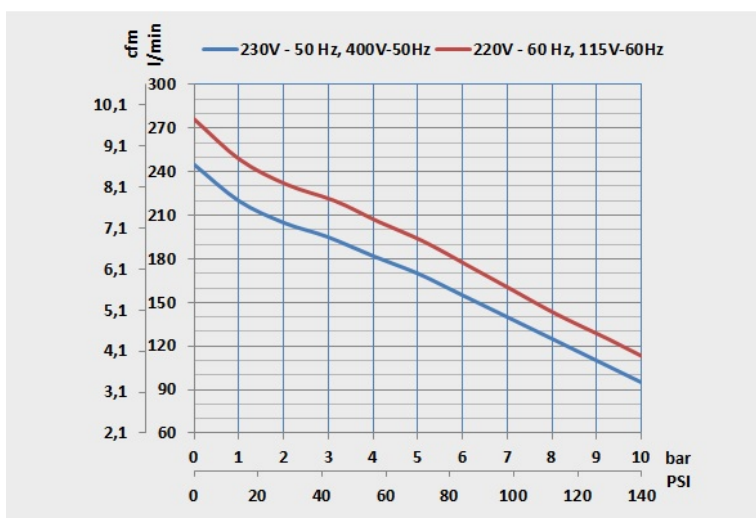

## Maintenance kits

Kit.	Description	Type	Hours
TA-FL-026-K-SCP	<p><b>PRIME 2-3 M filter cartridges kit</b></p> <p>Cartridge replacement kit for ordinary maintenance of the PURE AIR PRIME 2.3 M range.</p> <p>Comprehends:</p> <ul style="list-style-type: none"> <li>• 3 suction filters, code: TA-FL-058-C</li> <li>• 1 5 micron filter cartridge, code: TA-FL-001-C-SM</li> <li>• 1 0.01 micron filter cartridge, code: TA-FL-024-C2-SM</li> <li>• 2 o-rings for filters, cod. TA-FL-001-OR</li> <li>• 1 non-return valve pads, code: CO-PV-001</li> <li>• 1 5 micron filter cartridge (AF30), code: TA-FL-075-C</li> <li>• 1 0.01 micron filter cartridge (AFD30), code: TA-FL-076-C</li> <li>• 2 o-rings for filters, cod. TA-FL-075-OR</li> <li>• 1 AMG150 water separator filter cartridge, code: TA-FL-074-C</li> </ul> <p>Replace the filtration cartridges every 2 years.</p>	Ordinary	1.000h/1y
PR-GC-031-K-SCP	<p><b>PRIME 2</b></p> <p>Extraordinary maintenance spare parts kit PRIME 1.5 kW (PRIME 2) versions.</p> <p>comprehends:</p> <ul style="list-style-type: none"> <li>• 2 suction filters, code: TA-FL-058-C</li> <li>• 4 O-rings, code: MA-PC018</li> <li>• 2 piston rings , code: MA-MG00001</li> <li>• 2 piston guide rings, code: MA-MG00003</li> <li>• 2 expansion rings for the segments, code: MA-MG00024</li> <li>• 2 complete valve plates, code: MA-MG-02002P-C-P</li> <li>• 1 non-return valve pads, code: CO-PV-001</li> </ul> <p>The extraordinary maintenance intervals recommended for this compressor are 4000 hours/4 years.</p>	Extraordinary	4.000h/4y