




50/10 GENESI S










Cod. D-GE150-040

Oil free compressor, FAD at 5 bar 135 l/min, tank capacity 40 L, without air dryer.


 Maximum simultaneous use: 1

 Dental Chairs: 2

Medical Device CE 0425

											
Air Dryer	Air capacity	FAD F _{ad} @ 5bar	Motor	Voltage Hz	Receiver	Noise Level	RPM	Duty Cycle	max P.	Size	net Weight
No	m ³ /min l/min 0.182 182	m ³ /min l/min 0.135 135	kW 1.1	230V 50Hz	l 40	db(A) 65	1380	50 %	bar 7	mm 600x410x770	Kg 38

Air delivery performance

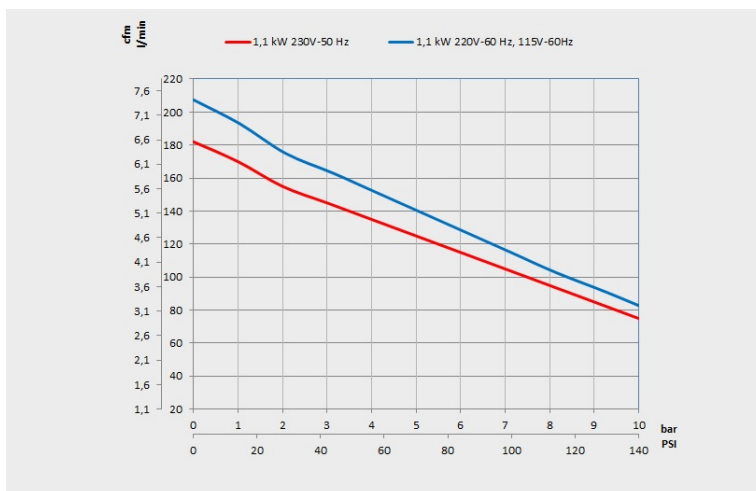

Maintenance kits

Kit.	Description	Type	Hours
TA-FL-001-K	GENESI S ordinary maintenance kit Ordinary maintenance kit for GENESI compressors, versions without dryer. The kit consists of the line filter cartridge (if present) and the pumping element suction filter. Annual replacement of kit components is recommended. Maintenance operations can be performed at a higher frequency if the environment in which the compressor works is particularly dirty.	Ordinary	1.000h/1y
MG-GC-011-K	Extra ordinary maintenance kit for GENESI S and M versions The GENESI extraordinary maintenance spare parts kit includes: a suction filter, two elastic rings, O-rings, a line filter cartridge and two complete valve plates. The recommended maintenance intervals for this compressor are 1500 hours / 2 years. Maintenance operations can be performed at a higher frequency if the environment in which the compressor works is particularly dirty.	Extraordinary	1.500h/2y