



## 50/10 GENESI M












**Cod. D-GE150-040-M**

Oil free compressor, FAD at 5 bar 120 l/min, tank capacity 40 L, with membrane air dryer and 0,01 micron filtration.


 Maximum simultaneous use: 1

 Dental Chairs: 2

Medical Device CE 0425

										
Air Dryer	Air capacity	FAD Fad @ 5bar	Motor	Voltage Hz	Receiver	Noise Level	Duty Cycle	max P.	Size	net Weight
Yes	m <sup>3</sup> /min l/min 0.182 182	m <sup>3</sup> /min l/min 0.135 135	kW 1.1	230V 50Hz	l 40	db(A) 65	50 %	bar 7	mm 710x410x770	Kg 42

## Air delivery performance

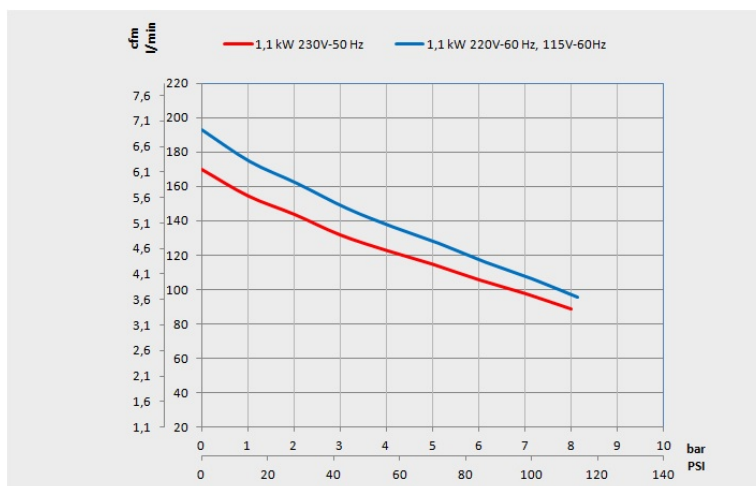

## Maintenance kits

Kit.	Description	Type	Hours
TA-FL-024-K-SCC	<p><b>GENESI M ordinary maintenance kit</b></p> <p>Ordinary maintenance kit for GENESI compressors with dryer. The kit consists of the cartridges of the dryer filtration battery (2 cartridges: 5 + 0.01 micron) and the relative O-rings, the pumping element suction filter and the non-return valve pad. Annual replacement of kit components is recommended. Maintenance operations can be performed at a higher frequency if the environment in which the compressor works is particularly dirty.</p>	Ordinary	1.000h/1y
MG-GC-011-K	<p><b>Extra ordinary maintenance kit for GENESI S and M versions</b></p> <p>The GENESI extraordinary maintenance spare parts kit includes: a suction filter, two elastic rings, O-rings, a line filter cartridge and two complete valve plates. The recommended maintenance intervals for this compressor are 1500 hours / 2 years. Maintenance operations can be performed at a higher frequency if the environment in which the compressor works is particularly dirty.</p>	Extraordinary	1.500h/2y